

Intelect® TRANSPORT ULTRASOUND

A New Dimension in Ultrasound


The Intelect Transport offers a new dimension in clinical ultrasound. The transportable ultrasound provides an innovative case design, with a logical control system and a large, easy to read graphical LCD display. User defined protocols are another great feature to customize any ultrasound application to your specific needs. Ergonomically designed ultrasound applicators are available in 2 cm², 5 cm² and 10 cm² as well as a new 1 cm² applicator for hand therapy. The new, repositionable base allows the Transport Ultrasound to be configured for desktop or wall mount use. The specially designed carrying bag (optional) can be used to protect and transport the unit and its accessories on the go.

Ultrasound

- Fully functional 1 and 3.3 MHz frequencies
- Ten new user-defined memory positions for User Protocols
- Pulsed and continuous therapy operation (10%, 20% 50% and 100%)
- Beams nonuniformity ratio (BNR) 5:1 maximum
- 110V or battery powered option
- New, ergonomically designed soundheads in variety of sizes (1 cm², 2 cm², 5 cm² and 10 cm²)
- Soundheads are interchangeable with Electronic Signature™ residing on all transducers
- Head warming and head contact control with visual and auditory feedback on all soundheads

Ordering Information

2782 – Intelect Transport
Ultrasound

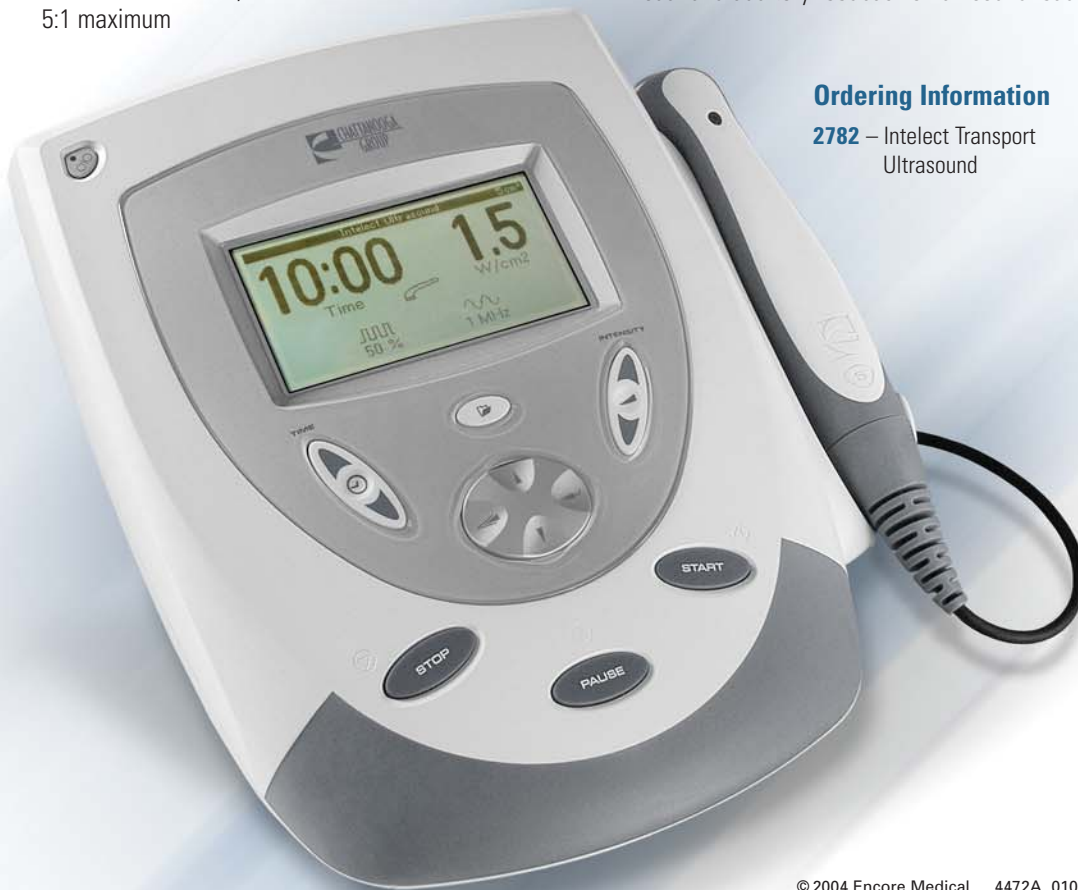

CHATTANOOGA
GROUP A DIVISION OF
encore
MEDICAL
ISO 13485 CERTIFIED